

Image 3



exclusive design to meet all your needs

Image 3 exclusive design to meet all your needs



Image 3 Standard



Image 3 Washable

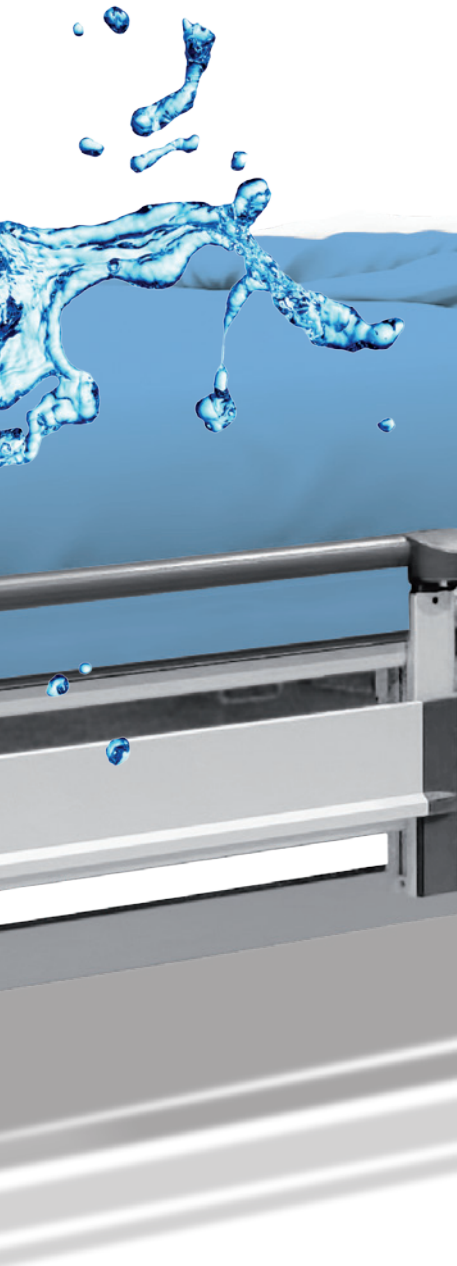


Image 3 XXL

Benefits of Image 3



extra
low position

28_{cm}



– Extra low position to prevent falls

– Support for independent mobilisation

– Unique modern design hand in hand with high comfort



Patient

– Sophisticated ergonomics for comfortable handling

– Automatic preset functions

– Easy bed hygiene



Staff

– Universal use – standard, bariatric and washable version

– Variable solution to meet the needs of individual clients

– High standard of infection control



Management

Meeting everyone's needs



Image 3 WASHABLE

The bed is available in a version suitable for cleaning and disinfection in automatic hospital washing machines.

- The bed construction undergoes a special technological modification that increases its service life and durability.

WASHABLE
FRIENDLY

Image 3 STANDARD

The modern design of the Image 3 bed will find its universal use in hospitals.

- By choosing the version, design and accessories, the bed can be adapted for a specific purpose.
- The attractive design of the bed helps to create a cozy and pleasant environment in health facilities.





Image 3 XXL

Bariatric version of Image 3 meets the growing demand for beds with higher load capacity.

- The bed has a wider mattress platform.
- The safe working load is increased to 320 kg.
- Its construction and accessories are safe even for a higher load.

WIDER BED

HIGHER LOAD

Variable design

HOSPITAL
DESIGN



EXCLUSIVE
DESIGN



UNIVERSAL
DESIGN



The design of the bed is timeless and suitable for a broad spectrum of needs in hospital care. Image 3 Standard, Washable and XXL fully respect the highest ergonomic principles in healthcare.



Due to its classic concept, the hospital design suits standard healthcare practice and in addition offers a version specially adapted for cleaning in automatic hospital washing machines.



The exclusive design is enhanced by a combination of aluminium elements and wooden decors. The elegant look has a comforting effect and creates a beautiful environment for both the patients and staff.



The simple design offers a compromise of the bed's unique practical features in an economic version.

Patient safety

Safety is a key parameter for the Image 3. All the bed's functions are designed so that the patient can not just feel, but really be as safe as possible.

Ergoframe®

The special mattress platform positioning significantly reduces the risk of pressure ulcers.



Hygiene


Thanks to the shape and material of its construction, the bed and its design elements are easy to clean.



BRAKE Alarm))


Brake alarm

A sound signal will notify the nursing staff of the need to brake the bed to ensure the client's maximum safety.

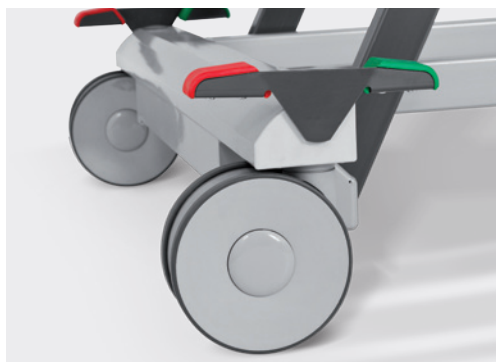

GO Button

Thanks to this patented safety concept, the bed's controls lock after three minutes of inactivity.

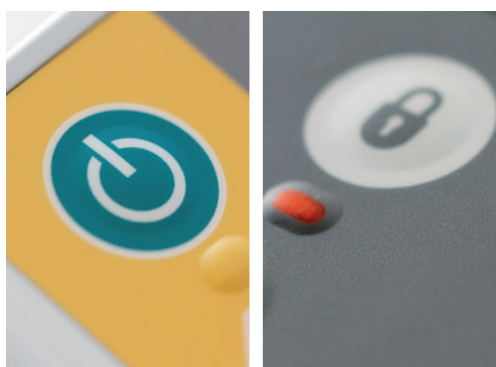

Functions lock

The electrical functions of the bed can be simply locked.

SafeFree[®] siderails

4 different adjustable heights of siderails. No required additional space beside the bed folding. Covered full length of bed sides



SafeSense

The SafeSense bed exit alarm warns the staff when the patient leaves the bed.



3-STOP strategy

EXAMINATION
POSITION



MOBILISATION
POSITION



EXTRA LOW
POSITION





Examination position

- The highest bed position provides the staff with ideal ergonomic conditions for work.
- Comfortable nursing procedures and care.
- The position provides complete access to the bed for cleaning and disinfecting.

80
cm



Mobilisation position

- Preset functions on the control panel for simple mobilisation.
- In combination with the siderail handle, the position simplifies and supports the mobilisation of the patient.
- The ideal setting for cleaning in an automatic hospital washing machine.

40
cm



Extra low position


- Separate button as a preset function on the control panel.
- Fall prevention.
- Safe position for sleeping.

28
cm



Practical details

- 1 5th castor
facilitates comfortable handling in all spaces
- 2 Bed extension
to three different levels according to the client's individual needs
- 3 Night light
for the patient's safe orientation
- 4 3D corner bumpers
protect the bed in a horizontal and vertical position
- 5 Mattress platform width
optional mattress platform width 86 or 90 cm with removable plastic parts
- 6 Linen shelf
may also be used to store the supervisor panel
- 7 Mechanical CPR
for emergency use

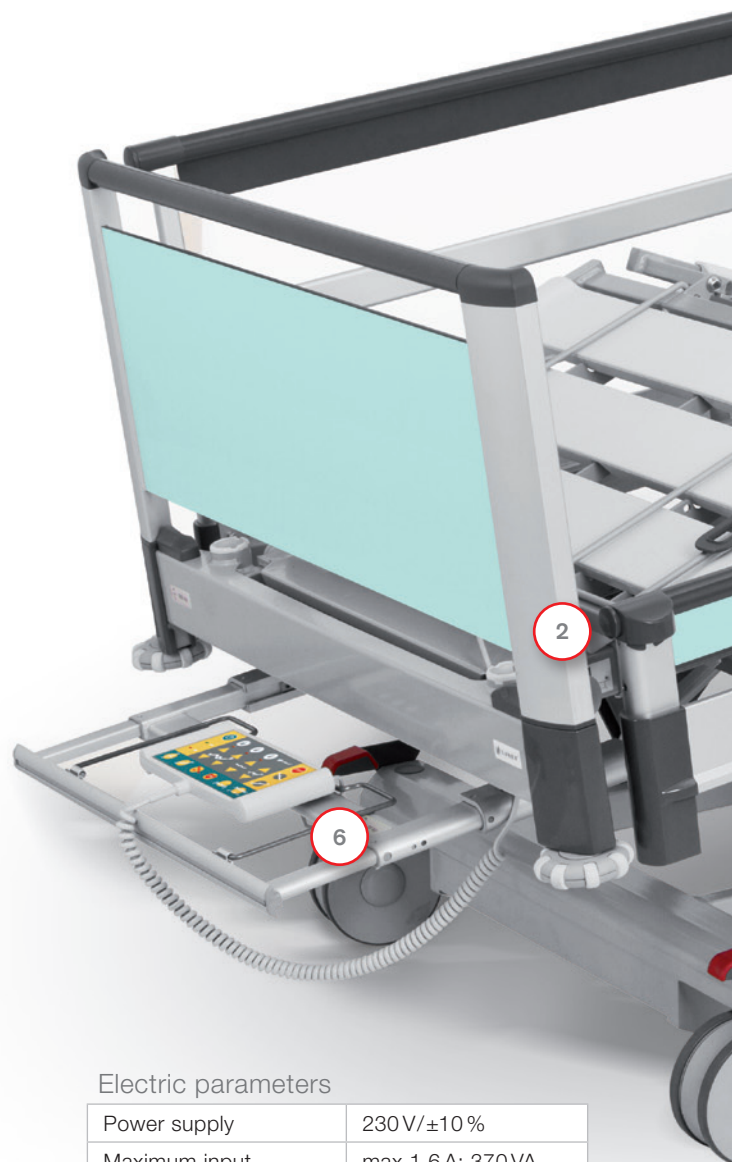

Electric parameters

Power supply	230 V/±10 %
Maximum input	max 1.6 A; 370 VA
Protection against water and dust	IPX4 (IPX6 for Image 3 Washable)

Certificates


LINET is certified according to European standards and the quality system:
ČSN EN 9001:2000 ČSN EN 13485:2003

The products are manufactured according to European standards for healthcare technology and hospital beds: ČSN 60601-1, ČSN 60601-2, ČSN 60601-4, ČSN 60601-2-52, ČSN EN ISO 14941
Image 3 Washable: AK-BWA:2009

Color concept

offers great variability: metal or plastic parts in silver or white color and a head ends with panels in decors of the customer's choice.



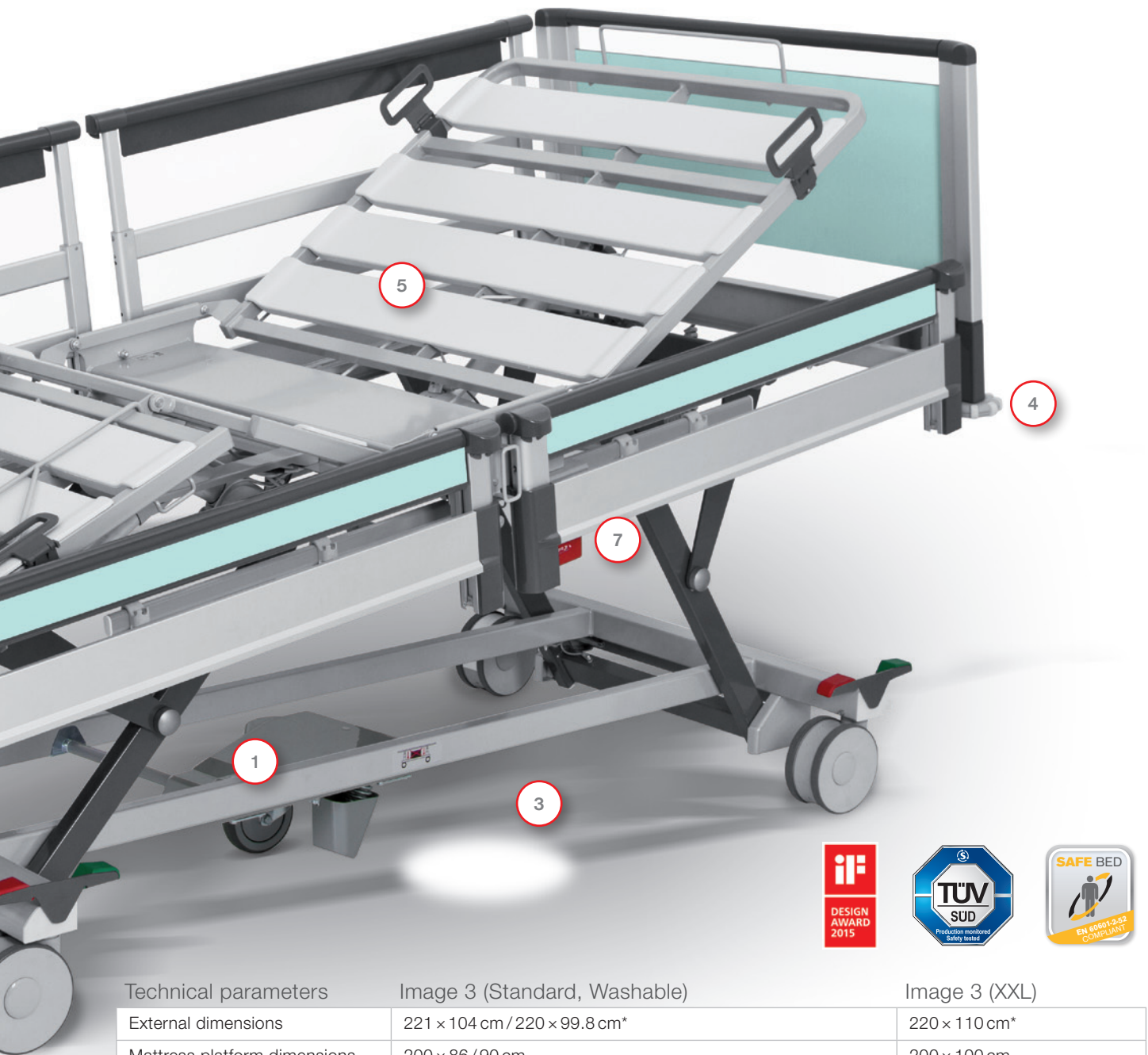

Technical parameters

Image 3 (Standard, Washable)

Image 3 (XXL)

External dimensions	221 × 104 cm / 220 × 99.8 cm*	220 × 110 cm*
Mattress platform dimensions	200 × 86 / 90 cm	200 × 100 cm
Clearance	13 / 15.5 cm	13 / 15.5 cm
Mattress platform height adjustment	28–80 cm	28–80 cm
Bed extension	10; 22; 30.5 cm	10; 22; 30.5 cm
Maximum mattress height	18 cm (split siderails); 17 cm (single collapsible siderails)	26 cm including the extender
Safe working load	265 kg (Standard) / 265 kg (Washable)	320 kg
Bed weight	160 kg *, **	160 kg *, **
Backrest tilt	0–70°	0–70°
Thighrest tilt	0–34°	0–34°
TR/ATR tilt	0–14°	0–14°

* by configuration; ** by accessories

Image 3

Patient

- Extra low position to prevent falls
- Support for independent mobilisation
- Unique modern design hand in hand with high comfort



- Sophisticated ergonomics for comfortable handling
- Automatic preset functions
- Easy bed hygiene



Staff

Management

- Universal use – standard, bariatric and washable version
- Variable solution to meet clients' needs
- High standard of infection control



 wissner-
bosserhoff

 LINET®

Members of LINET Group

LINET spol. s r. o.

Želevčice 5 | 274 01 Slaný | Czech Republic

tel.: +420 312 576 400 | fax: +420 312 522 668 | e-mail: info@linet.com | www.linet.com

