

# Eleganza 1 Extra SAFE

Maximum performance for specific departments



## Patient

- Safety
- Aid to reduce health risks



## Staff

- Easy bed hygiene



## Management

- Increased standard for specific unusual use




# High added value

The Eleganza 1 Extra SAFE is equipped with sophisticated functions which bring comfort, safety and aid in the treatment of patients and to healthcare staff.

## Ergoframe®

When adjusting the back and thigh rest the Ergoframe® system extends the mattress platform by 16 cm thereby minimising the pressure and friction in the patient's pelvic area. The risk of pressure ulcers in the exposed area is reduced.

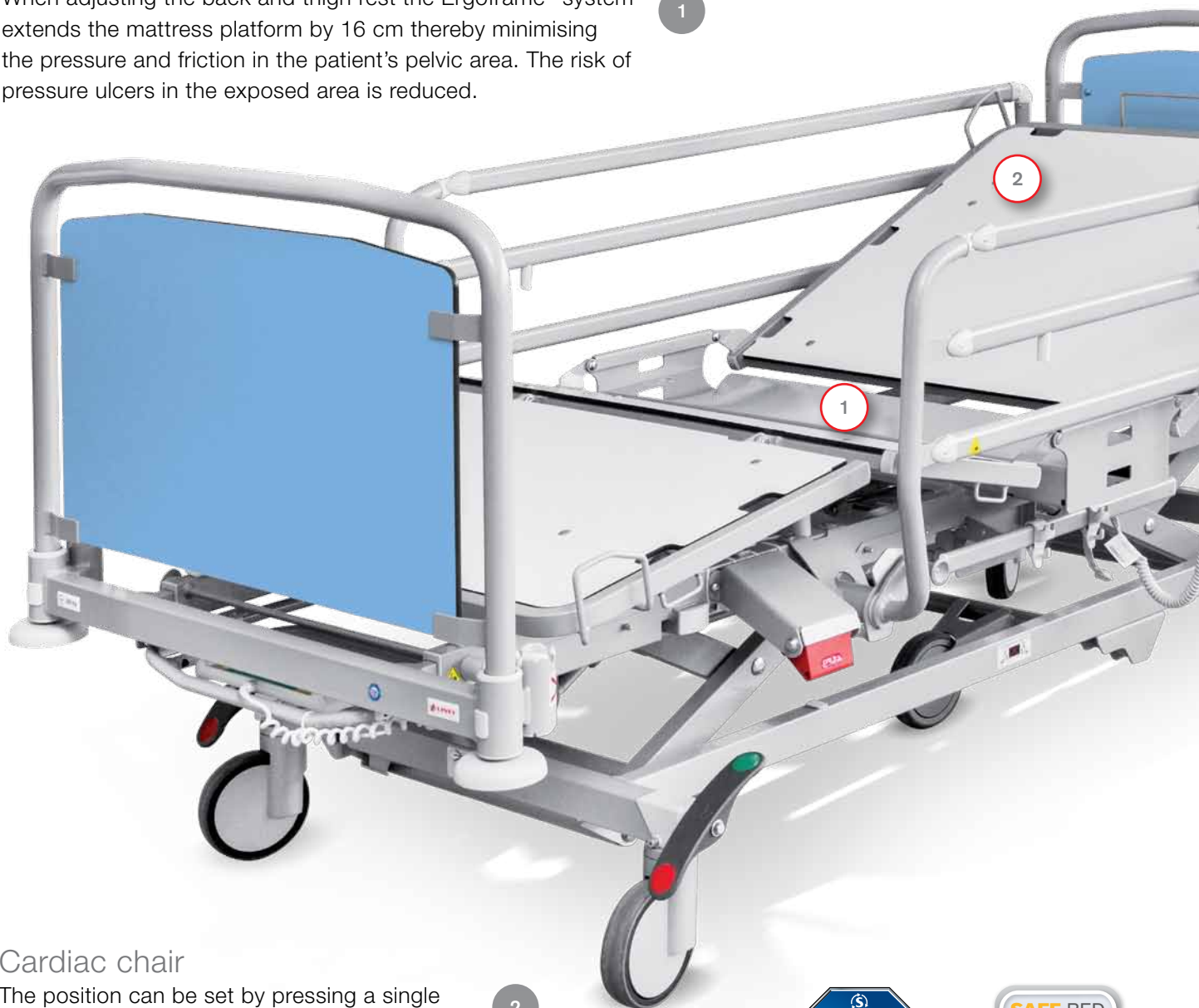

## Cardiac chair

The position can be set by pressing a single button. It helps to optimise the performance of the lungs and circulatory system. In addition, it is comfortable for the patient.



# Care and absolute safety


The Eleganza 1 Extra SAFE offers a high quality standard of hospital bed. Thanks to the modifications it is suitable for specific departments such as psychiatric, prison hospital, mental health etc.



Firmly installed bed ends and panelling of the ends. The mattress platform parts are fixed. All smaller parts of the bed, such as the hooks on the accessory bar, are secured.



Supervisor panel is built into the linen shelf.

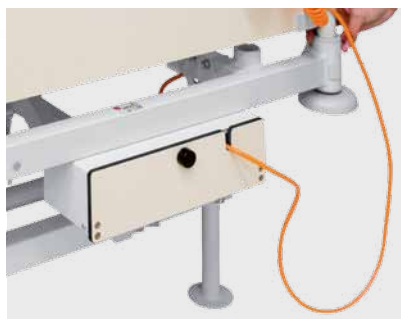

Power supply cable protection in a box with a lock. The cables are fed in a way so they do not endanger the patient.



SEGUFIX KEY is compatible, and the operator does not need any additional keys.



Bed in a version without castors can be fixed to the floor.

# Technical data

## Basic technical parameters

Safe working load (SWL)	185–250 kg
Bed weight (according to the accessories)	135 kg
External bed dimensions	218 × 99 cm
Internal bed dimensions (mattress size)	200 × 90 cm
Setting the mattress platform height	
without castors	36.7–73.7 cm
with castors 150 mm	39.5–77.0 cm
Clearance for servers	
without castors	16.0 cm
with castors 150 mm	18.8 cm
Anti-Trendelenburg	15°
Trendelenburg	15°
Painted metal bed parts	RAL 9002 – light grey

## Basic equipment of Eleganza 1 Extra SAFE \*

Fixed mattress platform covers
Fixed tube ends
Fixed end decors
Linen shelves
Inbuilt full supervisor panel
Plug and play connector of the patient's handset
Back-up battery
Supervisory panel with cardiac chair functions, TR and ATR by a single button
Protection of the supply cable in a locking box

\* Unlike the standard Eleganza 1 version



## Members of LINET Group

### LINET spol. s r. o.

Želevčice 5 | 274 01 Slaný | Czech Republic  
tel.: +420 312 576 400 | fax: +420 312 522 668 | e-mail: info@linet.com | www.linet.com

